



ADD-ON DRIVE
COMPATIBLE



A star is re-born!

The flagship model of the **Action_{NG}** family, is the stylishly, designed **Action³_{NG}**, based on the tried and trusted **Action³** chair, it offers improved functionality, comfort and stability. This medium active chair incorporates key design changes to enhance function with the use of lighter components to reduce weight providing an efficient rolling action. As with all models in the **Action_{NG}** family, all parts are cross compatible and interchangeable throughout the range whether you have an **Action_{NG} 2, 3 or 4** - the 'next generation' truly has arrived.



Yes, you can.®



Functionality

Invacare always seeks to innovate and the **Action³_{NG}** introduces a new style of cross bar with folding mechanism to improve stability and provides a smoother folding and unfolding action for easier handling. Plus it features a new design of footplate for greater comfort.

Durability

Modular and adaptable, the **Action³_{NG}** can accommodate progressive changes in customer needs with its wide range of options and accessories, which can be changed quickly and easily. With a proven track record for quality and durability, **Action³_{NG}** can be relied on for everyday use, day after day. This chair is also compatible with most of the **Alber** product range - please see the **Alber** brochures for further details.

Transportable

With a low transfer weight and foldable backrest, the **Action³_{NG}** is easy to fold and to dismantle, making it easily transportable in a car.

Features and Options



14.2kg weight

The total weight of 14.5kg addresses the need for improved handling, storage and transportation. The lowest transportation weight (without legrest, armrest and rear wheels) is 7.6 kg.



Brake and Brake extension

Fold the brake handle for easier side transfers.



Length compensating leg support

New elevating legrest offers optimal length compensation with the elevation stepped in approx. 3 – 5 degrees.



Footrest

The proven design of the swing in / swing out footrest has an easy-to-operate quick release mechanism. Adjust the foot support height in steps of 1 cm.

Features and Options



Tension adjustable backrest upholstery

Provides increased comfort and appropriate spinal support.



Amputee legrest

Height, angle and depth adjustable, the amputee legrest is also swingaway.



20 rear wheel positions

Driving performance is enhanced with a choice of 20 rear wheel positions as standard.



Handrims

Choice of plastic coated or anodized aluminium.



Reclining backrest

A gas spring backrest allows carers to easily adjust the back position in one continuous movement.

Features and Options



Hub brake

A hub brake is available for attendants. Push to lock brakes are included.



Transit chair

The Transit version offers seat heights of 46, 48.5 and 51 cm with 12" plastic rear wheels.



Buckle belt

Increased safety while seated.



New foot plate

A newly designed fixed footplate provides improved grip and comfort.



Front castor

Front castors available in 5", 6" and 8" pneumatic or solid tyre, various options of tyre width also available.



Head support

Adjustable on the backbrace to meet individual needs. Side and depth adjustment provides optimum positioning.

Features and Options



Armrest assortment

① Swing away detachable, ② swing away detachable height adjustable, ③, ④ and ⑤ removable height adjustable, ⑥ comfort detachable height adjustable and ⑦ hemi.



Anti-tip device

Provides extra safety on slopes or uneven ground.



Transparent tray

Useful for mealtimes and leisure activities such as reading.



One piece footplate

Optional one piece footplate adjustable in depth and angle.



Cane holder



Dual handrim propulsion


An improved dual handrim propulsion with QR mechanism makes the chair easily foldable.



Backrest types

On top of the fixed backrest (0° or 7°), **Action³NG** offers 4 different types of backrests to suit different needs :

- ① folding,
- ② angle adjustable (-6°/+24°),
- ③ racket reclining (0°/30°) and
- ④ gas spring reclining (0°/30°)

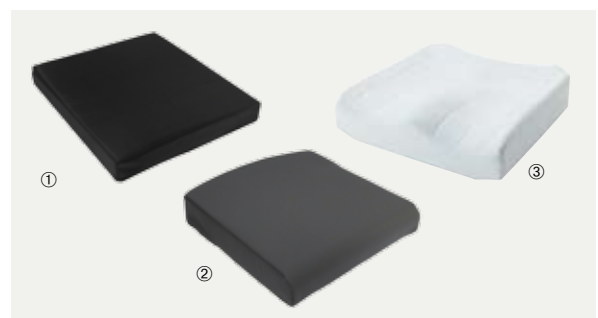

Cushions

① 5 cm foam cushion, ② comfort **Mistral2** cushion and ③ visco cushion.

Technical data

Invacare Action³ NG	38 / 40.5 / 43 / 45.5 / 48 / 50.5 cm	40 / 45 cm	38.5 * / 41 / 43.5 / 46 / 48.5 / 51 cm	40 / 43 / 51 cm	20 - 23 - 26 / 24 / 19 - 29.5 cm	35 - 48 / 34 - 47 cm	80° / 0 - 80°

Invacare Action³ NG	0° / 7° / 0° - 30° / -6° - +24°	Seat width + 20 cm	81.5 - 102 cm	100 cm	72 - 85 cm	28.5 cm	14.2 kg


Invacare Action³ NG	125 kg	Transport weight is 7.6 kg (without armrest, footrest and rear wheels)
		Crash tested ISO 7176-19**

* Available from January 2011
 ** Our wheelchairs comply above ISO standards and with ISO 10542, and have additional features providing increased levels of occupant security and safety whilst travelling in a motor vehicle. Wheelchair users should however transfer to the vehicle seat and use the vehicle installed restraint system whenever it is feasible.

Frame colours



Upholstery

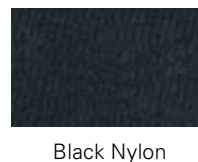

INVACARE® Limited
 Pencoed Technology Park - Pencoed -
 Bridgend CF35 5HZ
 Tel: +44 1 656 776 222
 Fax: +44 1656 776 220
 E-mail: uk@invacare.com
 www.invacare.co.uk
 Sales Order E-mail:
 ordersuk@invacare.com

INVACARE® Ireland Ltd
 Unit 5 - Seatown Business Campus
 Seatown Road - SWORDS - County
 Dublin
 Tel. +353 1 8107084
 Fax +353 1 8107085
 Email: Ireland@invacare.com
 www.invacare.ie

© 2010 Invacare® International Sarl
 All rights reserved.

All information quoted is believed to be
 correct at time of print.
 Invacare® reserves the right to alter
 product specifications without prior
 consultation.

Action3NG - UK - 11/2010
 1528014

