

INVACARE® TOP END® Eliminator OSR



Yes, you can.®

INVACARE® TOP END® Eliminator OSR




Edith Hunkeler - Switzerland

The *Invacare Top End Eliminator OSR* racing chair is designed for speed. The super-rigid roll-formed horizontal maintube, combined with a relaxed head tube angle for exceptional stability and aerodynamics, makes this model the choice of champion racers. If your sights are set on winning, the *Top End Eliminator OSR* will help you make it to the finish line, first!

The *Invacare Top End Eliminator OSR* is custom built to make the most of each athlete's potential by matching ability with the optimum position for maximum speed. The standard position is recommended for athletes who have limited balance and need more of a bucket or squeezed position. The athlete can sit in the traditional position with his/her feet on a footplate or use it in the kneeling position, which makes this model the ideal choice for beginners. The kneeling position is best for experienced athletes with good balance who only kneel in their chair or need an amputee seat.

Four different cage styles are offered to maximize each athlete's transfer skills while maintaining frame stiffness. The I-, V- and Open V-cage styles are offered on the standard OSR racing chair. The Open V-cage is the most popular because it allows more leg room and easier transfers. Only available on the kneel model, the U-cage requires perfect measuring in order to fit the racer correctly and optimize racing performance.

OPTIONS AND FEATURES

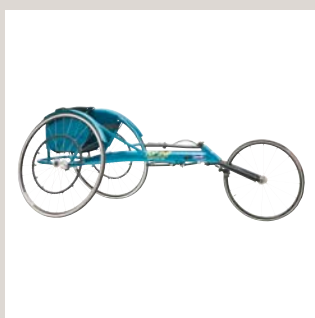


Top End Eliminator OSR featuring a relaxed head tube angle for exceptional stability.



Oversized rake comes as standard, Corima carbon front wheel available as an option.



The precision track compensator system ensures maximum speed on the road and around each turn of the track.

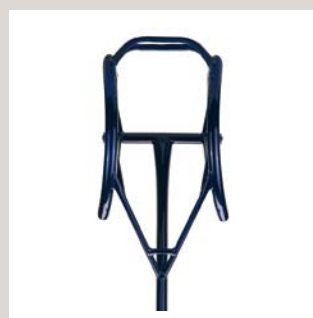

Open V-Cage shown with cut out camber tube option for more legroom.



TECHNICAL DATA



64" - 76"
(approx. 1626 - 1930 mm)



Max. 250 lbs
(approx. 113 kg)

Invacare Limited
Pencoed Technology Park - Pencoed
Bridgend CF35 5AQ - United Kingdom
Tel: +44 1 656 776 222
Fax: +44 1656 776 220
E-mail: uk@invacare.com
Sales Order E-mail:
ordersuk@invacare.com
www.invacare.co.uk

Invacare Ireland Ltd
Unit 5 - Seatown Business Campus
Seatown Road - SWORDS
County Dublin - Ireland
Tel. +353 1 8107084
Fax +353 1 8107085
Email: Ireland@invacare.com
www.invacare.ie

© 2012 Invacare International Sàrl
All rights reserved.

All information quoted is believed to be correct at time of print. Invacare reserves the right to alter product specifications without prior consultation.

TopEnd Eliminator - UK - 03/2012 - 1524412

Please refer to order form for all technical specifications.



Yes, you can.®